

Alcohol Testing Program - Instrument Processing Sheet

Agency: CLAY COUNTY SO

Instrument Serial Number: 80-001046

Date In: 8/27/2026

DI Completion Date: 8/28/26

☒ Ship

☐ P/U

☐ H/D

☐ CMI

☐ EE

Intake By: <u>WKP</u> Date: <u>8/27/2026</u>		Quality Checks By: <u>IAS</u> Date: <u>8/27/2026</u>		Flow Adjustment By: _____																
☒ Annual ☐ Dropped Off ☐ Registration ☐ Return from CMI / EE ☐ Training Instrument Visual Inspection ☒ Case ☒ Handle ☒ Keyboard ☒ Dry Gas Shelf ☒ Feet ☒ Breath Tube ☒ Ports ☒ Screws Tight Other Equipment/Accessories ☐ Power Cord ☐ Printer Cable ☒ Static Bag ☐ 12V DC Cable Notes:		☒ Breath Tube Screen ☒ Replace External O-Rings ☒ Instrument Set Up Verified ☒ R-Value: <u>135</u> ☒ Flow Verification (L/s) Flow Column #: <u>ATP102</u> 32 mm <u>.144</u> (.139-.169) 36 mm <u>.160</u> (.156-.190) 53 mm <u>.230</u> (.228-.278) 103 mm <u>.500</u> (.447-.547) ☒ Barometric Pressure Check Gauge ID #: <u>28427</u> Gauge: <u>1009</u> Instrument: <u>1011</u> ☒ Stability Checks <table border="1" style="width: 100%; border-collapse: collapse; margin-top: 5px;"> <tr> <th>Simulator</th> <th>Serial #</th> <th>Lot#/Exp</th> </tr> <tr> <td>0.050</td> <td>MP5088</td> <td>202603A 3/3/2028</td> </tr> <tr> <td>0.080</td> <td>MP5089</td> <td>202603B 3/3/2028</td> </tr> <tr> <td>0.200</td> <td>MP5090</td> <td>202603D 3/4/2028</td> </tr> <tr> <td>0.080 DGS</td> <td>N/A</td> <td>23725080A1 10/5/2027</td> </tr> </table>		Simulator	Serial #	Lot#/Exp	0.050	MP5088	202603A 3/3/2028	0.080	MP5089	202603B 3/3/2028	0.200	MP5090	202603D 3/4/2028	0.080 DGS	N/A	23725080A1 10/5/2027	Flow Column #: _____ ☐ 5L/min – 17mm ☐ 15L/min – 53mm ☐ 30L/min – 103mm ☐ R-Value: _____ ☐ Post Adjustment Verification (L/S) Flow Column #: _____ 32 mm _____ (.139-.169) 36 mm _____ (.156-.190) 53 mm _____ (.228-.278) 103 mm _____ (.447-.547)	
Simulator	Serial #	Lot#/Exp																		
0.050	MP5088	202603A 3/3/2028																		
0.080	MP5089	202603B 3/3/2028																		
0.200	MP5090	202603D 3/4/2028																		
0.080 DGS	N/A	23725080A1 10/5/2027																		
Optical Bench Adjustment By: _____		Department Inspection By: <u>KTS</u>																		
Barometric Pressure Gauge: _____ ID#: _____		Barometric Pressure ID#: <u>28427</u>																		
Simulator	Serial #	Lot #	Expiration	Gauge: <u>1010</u> Instrument: <u>1013</u>																
0.000		N/A	N/A	Mouth Alcohol Solution Lot #: <u>2025-D</u> Exp: <u>9/25/2027</u>																
0.040				Acetone Stock Solution Lot #: <u>2025-B</u> Exp: <u>9/22/2027</u>																
0.100				Simulator	Serial Number															
0.200				0.000	MP5086															
0.300				Interferent	MP5087															
0.080 DGS	N/A			0.050	MP5088															
☐ Post Optical Bench Adjustment Stability Checks				0.080	MP5089															
Simulator	Serial #	Lot #	Expiration	0.200	MP5090															
0.050				Attachments																
0.080				☒ Form 41 ☒ Stability Checks ☒ Calibration Certificate ☐ Optical Bench Adjustment ☒ Average Calculations ☐ Post-Stability Checks ☐ Flow Adjustment ☐ Form 40 ☐ Other:																
0.200																				
0.080 DGS	N/A																			
Gauge ID #: _____																				
Gauge: _____ Instrument: _____																				
Notes/Suggested Service:				☒ Instrument Complies with Chapter 11D-8, FAC ☐ Instrument Does Not Comply with Chapter 11D-8, FAC ☒ Return to/Place into Evidentiary Use ☐ Remain Out of Evidentiary Use ☒ Conduct an Agency Inspection Before Evidentiary Use																
				S. LeAndra Higginbotham Digitally signed by S. LeAndra Higginbotham Date: 2026.09.03 10:59:00 -04'00' Wen-Chi Pierson Digitally signed by Wen-Chi Pierson Date: 2026.09.08 09:49:24 -04'00'																
				Tech Review Admin Review																

Stability Checks

0.050 g/210L 0.047 to 0.053 g/210L	0.080 g/210L 0.077 to 0.083 g/210L	0.200 g/210L 0.194 to 0.206 g/210L	DGS 0.080 g/210L 0.077 to 0.083 g/210L ≤0.003 g/210L of Wet
Cup 210LW 80 Analyser - Manual Analyser Model 800 18/27/2008 Software: 8.000.27 Test 9/2/08 Time Air Blank 0.000 15:13 Control Test 0.049 15:14 Air Blank 0.000 15:14 Control Test 0.049 15:15 Air Blank 0.000 15:16 Control Test 0.049 15:16 Air Blank 0.000 15:17 Control Test Stats Average 0.047 Std Dev 0.000 Rel Std Dev(%) 0.000	Cup 210LW 80 Analyser - Manual Analyser Model 800 18/27/2008 Software: 8.000.27 Test 9/2/08 Time Air Blank 0.000 15:09 Control Test 0.077 15:10 Air Blank 0.000 15:11 Control Test 0.078 15:11 Air Blank 0.000 15:12 Control Test 0.078 15:12 Air Blank 0.000 15:13 Control Test Stats Average 0.0777 Std Dev 0.0016 Rel Std Dev(%) 0.7434	Cup 210LW 80 Analyser - Manual Analyser Model 800 18/27/2008 Software: 8.000.27 Test 9/2/08 Time Air Blank 0.000 15:41 Control Test 0.196 15:42 Air Blank 0.000 15:43 Control Test 0.196 15:43 Air Blank 0.000 15:44 Control Test 0.196 15:44 Air Blank 0.000 15:45 Control Test Stats Average 0.1960 Std Dev 0.0000 Rel Std Dev(%) 0.0000	Cup 210LW 80 Analyser - Manual Analyser Model 800 18/27/2008 Software: 8.000.27 Test 9/2/08 Time Air Blank 0.000 15:56 Control Test 0.075 15:58 Air Blank 0.000 15:59 Control Test 0.080 15:59 Air Blank 0.000 15:59 Control Test 0.079 15:59 Air Blank 0.000 16:00 Control Test Stats Average 0.0793 Std Dev 0.0006 Rel Std Dev(%) 0.7277
<i>Isabella Smyth</i> Operator's Signature	<i>Isabella Smyth</i> Operator's Signature	<i>Isabella Smyth</i> Operator's Signature	<i>Isabella Smyth</i> Operator's Signature

Florida Department of Law Enforcement Alcohol Testing Program

DEPARTMENT INSPECTION REPORT - INTOXILYZER 8000

Agency: CLAY COUNTY SO

Time of Inspection: 09:26

Date of Inspection: 08/28/2026

Serial Number: 80-001046

Software: 8100.27

Check or Test	YES	NO	Check or Test	YES	NO
Diagnostic Check (Pre-Inspection): OK	Yes		Date and/or Time Adjusted		No
Minimum Sample Volume Check: OK	Yes		Barometric Pressure Sensor Check: OK	Yes	
Alcohol Free Subject Test: 0.000	Yes		Mouth Alcohol Test: Slope Not Met	Yes	
Interferent Detect Test: Interferent Detect	Yes		Diagnostic Check (Post-Inspection): OK	Yes	

Alcohol Free Test (g/210L)	0.05g/210L Test (g/210L) Lot#:202603A Exp: 03/03/2028	0.08g/210L Test (g/210L) Lot#:202603B Exp: 03/03/2028	0.20g/210L Test (g/210L) Lot#:202603D Exp: 03/04/2028	0.08 g/210L Dry Gas Std Test (g/210L) Lot#:23825080A1 Exp: 10/05/2027
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.049	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079
0.000	0.050	0.078	0.196	0.079


Standard Deviations	0.0003	0.0000	0.0000	0.0000
---------------------	--------	--------	--------	--------

Average Standard Deviation of 0.05, 0.08 and 0.20 g/210L Tests: 0.0000 Number of Simulators Used: 5

Remarks:

The above instrument complies (☒) does not comply (☐) with Chapter 11D-8, FAC.

I certify that I performed this inspection in accordance with the provisions of Chapter 11D-8, FAC.

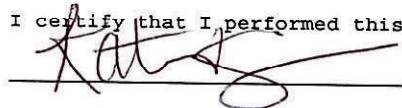

KATIE T SPEARIN
Signature and Printed Name

08/28/2026
Date



Calibration Certificate

Florida Department of Law Enforcement
Alcohol Testing Program
2331 Phillips Road
Tallahassee, FL 32308

This is to certify the calibration of Intoxilyzer 8000 serial number 80-001046 , manufactured by CMI, Inc. was calibrated in accordance with FDLE/ATP Form 36 - Department Inspection Procedures - Intoxilyzer 8000.

Serial Number:	<u>80-001046</u>	Average of $y \pm$ Uncertainty*
Owning Agency:	<u>CLAY COUNTY SO</u>	0.050 g/ 210 L 0.004
Calibration Date:	<u>08/28/2026</u>	0.078 g/ 210 L 0.004
Calibration Time:	<u>09:26</u>	0.196 g/ 210 L 0.008
		0.079 g/ 210 L Dry Gas Control 0.005

All results are reported in g/ 210 L.

Bias is limited by calibration acceptance criteria. All calibration results must be within ± 0.005 or 5%, whichever is greater, of the target alcohol concentration.

*Uncertainty is based on fleet-wide data and is expressed to a 99.73% level of confidence ($k=3$).

The instrument results before and after any adjustment are found in the associated pre and post stability checks.

TRACEABILITY INFORMATION

This instrument was calibrated using solutions prepared by Alcohol Countermeasure Systems, Inc. (ACS). ACS prepared and certified these CRMs in accordance with ISO 17034 and ISO/ IEC 17025 Standards.

Simulator temperatures are traceable to NIST. Simulator temperatures are checked with NIST traceable digital thermometers calibrated by Precision Metrology in accordance with ISO/ IEC 17025 standards.

Dry gas control measurements are traceable to NIST through the use of CRMs supplied by an accredited CRM supplier. The supplier of dry gas standard controls prepared and certified the CRMs in accordance with ISO Guide 34 and ISO/ IEC 17025 standards.

This document shall not be reproduced except in full, without written approval of the Florida Department of Law Enforcement Alcohol Testing Program.

Digitally signed by Kaitlyn
Spearin
Date: 2026.08.31 11:18:17
-04'00'

Kaitlyn Spearin

08/28/2026

Date

KATIE T SPEARIN,

Department Inspector

FDLE/ATP Form 69 August 2026
Issuing Authority: Alcohol Testing Program

Service • Integrity • Respect • Quality